™
CASADA
healthcare

< relaxmyway

BODYSHAPE
— Limited Edition —

< relaxmyway

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< relaxmyway

Since the company was founded in the year 2000, the name CASADA has stood for high grade products in the area of wellness and fitness.

Today, CASADA – products are available worldwide. What distinguishes CASADA products is their continuous development, unmistakable design and their absolute focus on improving quality of life as well as their reliable and excellent quality.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal BodyShape Limited Edition.

Table of Contents

12	Safety instructions
12	Scope of delivery
13	Technical Data
13	Assembly / Start-Up
14	Handling
14	Specials
15	Remote control
16	Troubleshooting
16	Storage / Maintenance
17	EU Declaration of Conformity
17	Warranty conditions

Safety Instructions

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- This device does not require any maintenance. Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- To avoid the danger of short-circuits, disconnect the power cord when not used for long periods of time.
- Do not use any damaged cables, plugs or loose sockets.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Misuse or incorrect use excludes any liability for damage.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- Do not use this device while driving.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.
- This equipment may be used by children of at least 8 years of age and by persons with reduced physical, mental abilities and a person with no experience or knowledge of the equipment, if supervision or instructions regarding the use of the equipment have been ensured in a safe manner so that the threats were understandable. Children should not play with the equipment. Unattended children should not perform cleaning and maintenance of equipment.

Scope of Delivery

- Power supply unit 24 Volt / 12 Volt
- Manual
- Massage mattress

Technical Data BodyShape LE

Measurements:	176 x 61 cm
Weight:	2,1 kg
Voltage:	input: AC 100-240V ~ 50/60Hz, 800 mA output: DC 12V, 2500 mA
Nominal power:	max. 30 watt
Automatic runtime:	10 minutes
Certificates:	

For your own safety we equip the BodyShape massage mattress with a protection against overheating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk and after a runtime of 10 Minutes, the BodyShape Limited Edition will switch off automatically and will only allow switching on again after an appropriate downtime.

Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 30 Minutes.

Assembly / Start-Up

1. Open the carton, take out the mains adapter and place it in reach.
2. Connect the provided adapter to the massage mattress and then connect the adapter to the socket.
3. To switch on the massage mattress, press the ON/OFF button on the remote control. The LED lamp lights up red and the device starts directly in the automatic mode M1 in medium intensity.

Handling

AUTOMATIC PROGRAMMES

The automatic mode massages in 5 different itineraries. Each program (M1 – M5) runs for a duration of 10 minutes and is selected via the „Mode“ button. During the automatic programmes you can choose between the 3 intensities anytime and switch the heating on and off. After finishing the selected program, the device automatically turns off.

MANUAL HANDLING

Via the remote control you can also select only one of the five different vibration fields. Doing this, the selection of intensity and warmth is still possible anytime. After the recommended operating time of 30 minutes, the device turns off independently.

Specials

VIBRATION

The vibration sets entire areas of the body in even oscillations which have far-reaching effects. Stimulating for the circulation and relaxing for the muscular system.

HEATING

A deep penetrating and heat-radiating effect in which one can perfectly relax. This function is particularly beneficial for muscle tension due to the deep reaching effect.

AUTOMATIC

5 specifically developed automatic programmes for the whole body, which offer a big range of different massage techniques. Additionally you can also turn on the heating function, if you'd like.

CAR ADAPTER

The handy car adapter makes your massage device portable.

The massage mattress comes with 10 intelligent build-in vibration fields, split into 5 zones, as well as 4 heating fields, which are spread in two zones each.

Remote Control

MESSAGE AREAS


NECK	2 vibration fields for a neck massage
UPPER BACK	2 vibration fields for a massage in the upper back area
LOWER BACK	2 vibration fields for a massage in the lower back area
LEG	2 vibration fields for a massage in the leg area
FEET	2 vibration fields for a feet massage

AUTOMATIC PROGRAMMES

M1	Vibration fields massage simultaneously in an up- and down going rhythm
M2	Long vibration intervals in alteration UPPER BACK – LOWER BACK – LEG – FEET – NECK
M3	Short intervals of individual and combined vibration fields NECK & UPPER BACK – LEG & FEET – LOWER BACK, afterwards all fields mutual in a fast rhythm
M4	Fast, short vibration intervals in alteration UPPER BACK & LOWER BACK – NECK & FEET
M5	Vibration fields massage simultaneously in an up- and down going wavy rhythm

FUNCTIONS

ON/OFF	With this button you can either switch the massage mattress ON or OFF
HEAT	ON and OFF switch of the heating fields
INTENSITY	Adjustable intensity in 3 levels
MODE	Selection of the particular automatic programmes

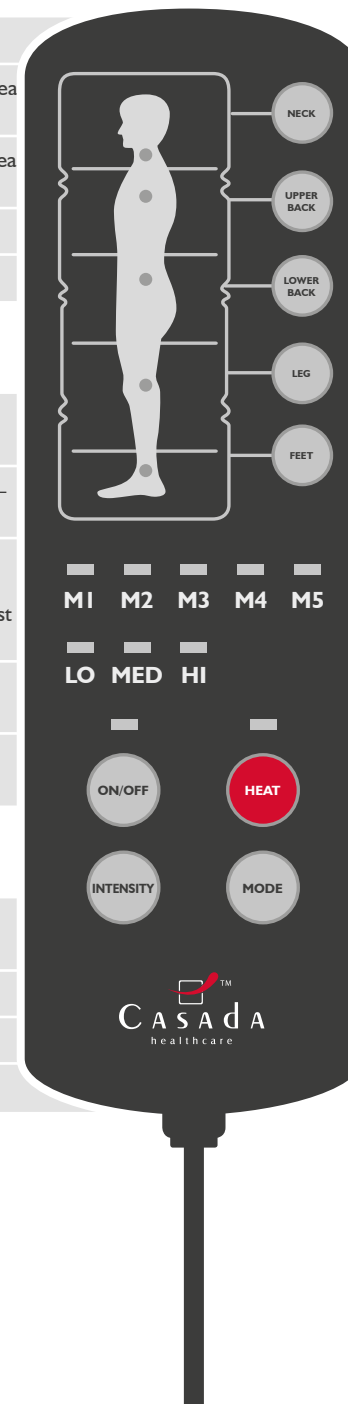Troubleshooting

Fault	Remedial action
The device is connected to the mains, however neither the device nor the remote control operation work.	<ul style="list-style-type: none">• Check whether all connections are correctly paired. Ensure that the power supply is fully connected.• The device can overheat if used permanently. Therefore, turn it off, disconnect it from the mains power supply and let it cool down completely.

Storage / Maintenance

With regular maintenance you can increase the service life of the BodyShape Limited Edition and ensure safe operation at all times.

1. Recommended storage: Dry, free of dust.
2. Keep liquids or open flames away from the device.
3. Clean the device with an imitation leather care and with a dry soft cloth.
4. Disconnect the device from the main power supply when it is not in use.

EU Declaration of Conformity

Compliant with European regulatory standards EMC 2014/30/EU, LVD 2014/35/EU.

Warranty Conditions

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer. The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medical treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.
Copyright © 2018 Casada International GmbH.
All rights reserved.

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Copyright © 2018 CASADA International GmbH.